

# Better Service Better Performance

## Technical Parameter

Type of Cassette: Standard

Capacity of Loading : 50,75,100 cassettes as required.

Output capacity of collection: Standard: 10 cassettes

Optional tray: 64 cassettes

Hooper: 1pc

Two type hopper option: Plastic for sleeve package

Metal for taped package

Print Speed: 3 seconds/cassette (20 cassettes/minute more)

Print Resolution:2500dpi

Date interface:USB

Data Acquisition: Barcode Reader

Voltage:110-240VAC

Power Draw:120va

Dimension(WxDxH):220\*400\*365mm

Weight:19kg



For Your Better Experience

## Laser Cassette Printer

PrintSmart LCPlus-S



TRIA-LINK

**Changzhou Trialink Co.ltd**

Add:1408, Building 3, Manhattan Plaza, Xinbei District, Changzhou City, Jiangsu, China

Website: <http://www.tria-link.com> Email : [info@tria-link.com](mailto:info@tria-link.com)

**Trialink Co.ltd**

Add: Unit G25 Waterfront Studios, 1 Dock Road, London,UK



## ADVANCED TECHNOLOGY


- Non Contact laser to cassette marking technology  
Print identification information on the cassette to create permanent directly
- Displace ink cartridges or ribbons, consumable free
- The operation cost is greatly reduced when compared with the traditional printers.
- Minimal Maintenance and Easily to do.

Various printing characters:

texts, symbols, graphics, logos, 1D&2D codes, barcodes for reliable cassette tracking.

Work on-demand and batch printing

## PrintSmart LCPlus-S

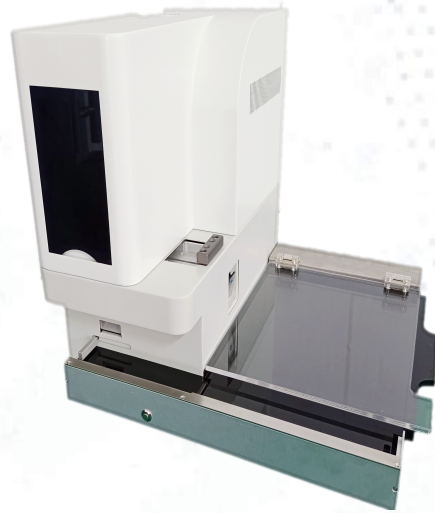

Smart delivery system protect cassette and safety transportation.

Reduce human errors and flexibility working.

Clear print quality and permanent labels for histology long life storage.

Reliable identification to enable the best tracking capacity in your laboratory.

## DURABILITY AND RELIABLE

## FRIENDLY LINK TO WORKFLOW

Compatible with most of color cassettes in market.

Compatible with LIS and HIS

Compatible with barcode scanner.



## STABLE PERFORMANCE

Easy Operation: Touch screen option

Safe Operation: Carbon Filter and Unique design for air purification

Friendly Operation: Smart software to flexible template edition

Language: English or Multiple languages available


## KEY PARAMETERS

More option of Hoopers: Metal with 50 cassettes, 75 cassettes, 100 cassettes or plastic sleeves.

Standard collection tray could face to left or right side as requirement.

Optional Collection Tray for 64 cassettes collection .

Smaller Size to suitable any place.

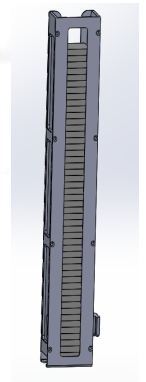

TRIA-LINK