



DAKEWE

CS600

Glass Coverslipper

Enjoy the smooth coverslipping process



Enjoy walkaway process with massive capacity

500-coverslip large magazine continuously loadable and unloading station holding up to 300 slides free technicians from repetitively loading and unloading and allow for an uninterrupted workflow.




Easy to use

A sleek and intuitive interface, delivering a seamless and effortless user experience.



Slide traceability

Equipped with an integrated barcode scanner, CS600 connects to Dakewe stainers to deliver complete process monitoring - tracking protocols, reagents, and operations while enabling seamless LIS integration.

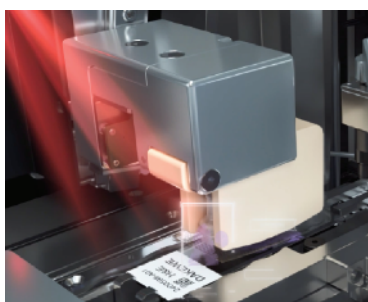

Care for the users from the air they breathe

CS600 coverslipper is equipped with a unique internal airflow circulation system, which ensures effective fume filtration and guarantees optimal user safety.



Hassel-free with all statuses in control

Precision level sensors continuously monitor mounting medium and cleaning reagent volumes, preventing process interruptions and eliminating the urgency of last-minute reagent replenishment.

Unsupervised worry-free



coverslipping



Adaptability



Seamlessly handles 50, 55, and 60 mm coverslips with adaptable magazines while offering intuitive mounting medium control for all applications, compatible with most commercially available mounting media.

Fast drying




Heated and well ventilated uploading station accelerates slide delivery.

Versatility



The coverslipping procedure is highly flexible, featuring both dry and wet coverslipping options to accommodate diverse customer needs. The mounting medium volume can be adjusted in real time to suit different specimen types.

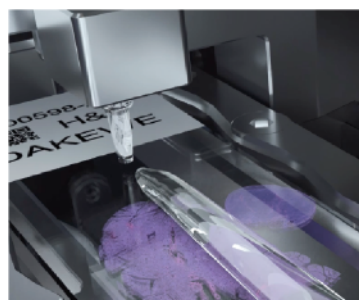

Strict coverslips quality control



Broken and adhered coverslips are detected by the dedicated sensor and removed to the defect tray, ensuring a seamless coverslipping process without human intervention.

Optimal coverslipping outcome



Automatic nozzle soaking and regular piping cleaning reduce the risk of blockage and bubbles, and further offer a high level of stability.

Specifications

Dimensions, (W × D × H)	700 × 650 × 665 mm
Net weight	Approx. 89 kg
Power consumption	250 VA
Supply voltage/frequency	100 - 240 VAC, 50/60 Hz
Slide dimension	(24.8 - 26.2) × (74.8 - 76.2) mm
Coverslip length	50, 55, 60 mm
Coverslip storage capacity	Up to 500 coverslips in 5 magazines
Unloading station storage	Up to 300 slides
Oven temperature	50 °C
Workstation	Can be used stand-alone or integrated with DP360 without extra connector

DAKEWE

Dakewe medical
Shenzhen, China
www.dakewemedical.com
Ver. 1.0, Rev. Mar 2025

