



DAKEWE

SurePrint S10/S10Pro
slide printer

Pocket-sized genius



SurePrint S10 slide printer

Incredibility in the tiny body

S stands for... speed, saving, and slim.

Maximized compactness

Got a space for a smart phone? Then there is a place for S10

With 115×145 mm minimum distance between feet, S10 can effortlessly save scarce space in laboratories

Minimized consumable cost and maintenance

UV laser marking eliminates the demand for printing consumables such as ribbons, print heads, and ink cartridges, and hence the relevant maintenance

Mellow user experiences

Optimized user interface in simple/advance mode suitable for both geeks and pragmatists

Operate directly on the touch screen as a standalone printer, or connect it to a computer or various LIS

Multi-to-one control:
Integrate S10 to your intranet to receive print requests from multiple computers for optimized equipment utilization

Relative size compared to Dakewe MT1 microtome



The pocket-sized genius

Dakewe SurePrint S10/S10Pro slide printer



Easy loading


150 slides capacity

Insert-n-print slides with tissues using dedicated on-demand loading hopper



Built-in scanner facilitates lean workflow

On-demand output can be upgraded to batch output with extension (only S10Pro)

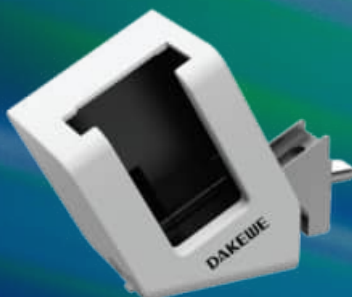

Specifications

	S10	S10Pro
Print technology	Non-contact UV laser marking technology	
Print content	Numbers, characters in multiple languages, symbols, 1D/2D barcodes	
Print speed	Up to 2 sec/slide, depending on print content	
Print resolution	2500 dpi	
Slide dimensions	(74.8 - 76.0) * (24.8 - 26.0) * (0.8 - 1.2) mm	
Load capacity	Up to 150 slides in 1 hopper	
Output	Up to 10 slides	Up to 50 slides
Built-in scanner	Yes	Yes
Air purification	Built-in air purification system	
System interface	LIS, HIS, PIC, PACS	
Operating systems	Windows 7, Windows 10, Windows 11, Linux	
Power	100 - 240 VAC, 50/60 Hz, 170 VA	
Dimensions	175 × 335 × 340 mm	175 × 440 × 340 mm
Weight	Approx. 13.5 kg	Approx. 15 kg

DAKEWE

Dakewe medical
Shenzhen, China
www.dakewemedical.com
Ver 1.0, Rev. July 2024

